coline[®]

Bathroom scale

Personvåg • Personvekt

Henkilövaaka • Personenwaage

Art. No.:

34-2163, 34-2220

Model. No:

BG-1003

English	2
Svenska	4
Norsk	6
Suomi	8
Deutsch	10

ENGLISH

SVENSKA

NORSK

SUOMI

DEUTSCH

Bathroom scale

Art.no. 34-2163, 34-2220 Model BG-1003

Please read the entire instruction manual before use and save it for future reference. We reserve the right for any errors in text or images and for making any necessary technical changes to this document. If you should have any questions concerning technical problems please contact our Customer Services.

Safety

- This product is a precision instrument and should be handled with care and protected from blows and rough treatment.
- Place the scale at the same location every time you weigh yourself.
- For accurate readings, the scale should be placed on a hard, level surface. An uneven floor beneath the scale can cause the readings to fluctuate.
- Weigh yourself without clothing at the same time of day, preferably in the morning before breakfast. Place your feet in the same place every time you weigh yourself.
- The displayed readings are rounded up/down to the nearest 100 g.
- The scale can be cleaned using a damp cloth, if and when required. Make sure that no moisture gets inside the product. Never use solvents or strong detergents.
- Keep the scale away from water and splashes, the electronic components can be damaged.
- Do not store the scale standing on its edge as the display will be on all the time and quickly drain the batteries.
- Remove the batteries from the scale if it is not to be used for an extended period.
- The scale can be affected by changes in temperature or by being moved. Inaccurate readings may occur if the scale is too cold when used. Allow the scale to stand and acclimatise for approx. 24 hours before using it for the first time.
- This product must not be used by anyone (including children) suffering from physical or mental impairment. This product should not be used by anyone who does not have the proper experience or knowledge of the product unless they have been instructed in its use by someone who takes responsibility for their safety.
- Never let children play with the product.

Batteries

The battery compartment on the bottom holds four AAA/LR03 batteries.

1. Open the battery compartment and insert the batteries. Refer to the markings in the battery compartment.
2. Refit the battery cover.


Note: “Lo” will appear on the display when the batteries need changing.

Operation

Warning: Bear in mind that the glass surface might become slippery when wet. Make sure that the surface is dry before stepping onto the scale.

1. Press [UNIT] (see fig.) to select the unit of measurement. Readings can be given in either kilograms (kg), pounds (lb) or stones (st).
2. The scale should be placed on an even level surface.
3. Carefully mount the scale with your weight evenly distributed between both feet. Remain stationary until the display shows a steady reading. The display flashes once the reading has been registered.
4. The scale automatically shuts off as soon as you step off it.

Note: If the scale is overloaded “Err” will appear on the display.

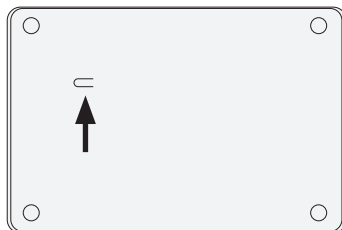Care and maintenance

Clean the product by wiping it with a soft, damp cloth. Never use solvents or abrasive cleaning agents.

Disposal

This product should be disposed of in accordance with local regulations. If you are unsure how to proceed, contact your local authority.

Specifications

Max weight	150 kg, 330 pounds, 24 st
Accuracy	100 g, 0.2 pounds
Batteries	4 × AAA/LR03 (included)

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