TENS/EMS device

TENS/EMS-apparat


TENS/EMS-laite

Quick guide

Snabbguide

Hurtigguide

Pikaopas

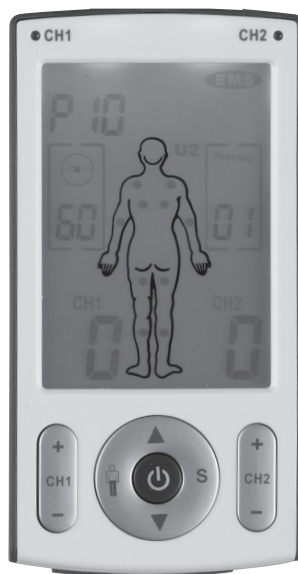Art.no

34-4103

Model

LT3011

Ver. 201010

English

Svenska

Norsk

Suomi

clas ohlson

TENS/EMS device

Art. No 34-4103 Model LT3011

Quick guide

Please read the entire instruction manual before use and save it for future reference.

We reserve the right for any errors in text or images and for making any necessary technical changes to this document. If you should have any questions concerning technical problems please contact our Customer Services.

1. Safety

- This product is intended for home use.
- Read and follow the user advice in these instructions. Bear in mind the restrictions and risks that accompany this device. Read the cautionary texts and symbols on this device.
- This product may not be used in the following circumstances:
 - Confirmed or suspected heart defects.
 - Some persons may experience skin irritation or hypersensitivity from electric stimulation or from the use of gels.
 - If a rash appears or pain persists, suspend treatment and consult a doctor.
- Placement of electrodes and settings should be based on advice from a doctor or a person who is used to administering treatment involving similar devices.
- Local skin irritation may occur from prolonged use of electrodes in the same position.
- Make sure the electrodes have good contact with the skin.
- If the stimulation level is uncomfortable and too high for you the impulses must be reduced to a more agreeable level. Suspend treatment and consult a doctor if you are unsure.
- This product should not be used by persons regularly out in traffic or persons who monitor machines and processes.
- Do not use this product on premises where aerosols or oxygen are stored or used.
- Do not use this product in the vicinity of explosive, highly flammable substances or gases.
- Do not use this product at the same time with other products that send electric impulses into the body.
- Do not connect electrodes to any other device or connect any other device, i.e., earphones and the like, to the TENS/EMS device.
- Do not use sharp objects to press buttons.
- Make sure that electrodes and cables are undamaged before using.
- The TENS/EMS device may only be used in conjunction with the accompanying accessories.
- Treatment may not be applied above the neck, i.e., to the face, on the throat, or anywhere on the head.

This product must not be used in the following circumstances

- Pregnant women should consult their doctors before using.
- This product should not be used if carcinogenic changes are found in the area treated.
- Treatment should not be applied to swollen, infected, inflamed or damaged skin.
- This product should not be used on persons who have a pacemaker or defibrillator.
- This product should not be used on persons who have diabetes, abnormally high blood pressure, or on persons who have epilepsy.
- This product should not be used on persons who have serious circulation problems in the body's lower parts, or hernia in the groin or abdominal area.
- In cases of known myocardial inflammation.
- Do not use the device for symptoms that have not been diagnosed.

Caution! Do not use the device if you are unsure how it works. Consult with a doctor or a person who is used to administering treatment involving similar devices if you feel unsure.

2. Inserting batteries

1. Open battery cover (on back) by applying pressure to it and at the same time push cover in the direction of the arrow.
2. Insert 4x AAA/LR03-batteries into battery compartment as indicated by markings.


3. Attaching cables to electrodes

Insert cable contacts in electrode contacts. Squeeze the contacts together so that no metal is visible at any connection.

4. Attaching electrode cables to the TENS/EMS device

1. First make sure that the device is switched off.
2. Grasp an insulated part of the contact and connect it to one of the outlets (15) and (16) on the top side.

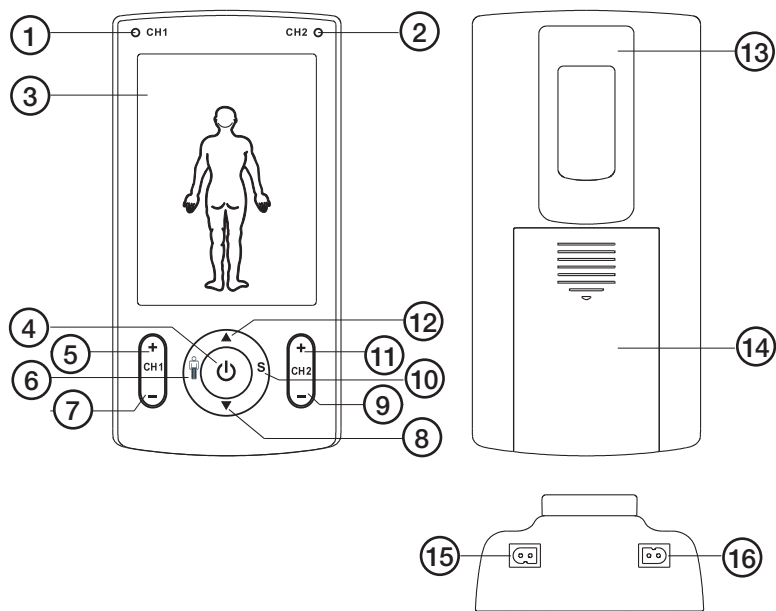5. Description

1. **CH1** lights up when channel 1 is activated.
2. **CH2** lights up when channel 2 is activated.
3. **LCD display**
4. **[⏻] Switch**. Press to switch on, press and hold down to switch off.
5. **[CH 1+]** – Increase intensity for channel 1. From **0** (switched off) to **60** (max).
6. **[U]** – Select user. Select **U1** (user 1) or **U2** (user2).
7. **[CH 1-]** – Switch off button lock. Reduce intensity for channel 1. From **60** (max) to **0** (switched off).
8. **[▼]** – Programme selection, back. **EMS** P20 – P01, **TENS** P20 – P01, **U** (own settings), U10 – U01, **MASSAGE** P10 – P01.
9. **[CH 2-]** – Switch off button lock. Lowers intensity for channel 2. From **60** (max) to **0** (switched off).
10. **[S]** – Press and hold for approx 3 seconds to set user parameters.
11. **[CH 2+]** – Increases intensity for channel 2. From **0** (switched off) to **60** (max).
12. **[▲]** – Programme selection forward.
13. Belt clasps.
14. Battery cover.
15. Connection for channel 1.
16. Connection for channel 2.

6. Buttons and functions

	<ol style="list-style-type: none"> 1. Switch for switching on the device and to start treatment. 2. Press to stop treatment. 3. Press and hold down button 3 seconds to switch off.
	Select treatment programme forward. EMS (P01 – P20) – TENS (P01 – P20) – U (own settings U01-U10) – MASSAGE (P01 – P10).
	Select user 1 or user 2.
[S]	Press button for pause and to resume treatment
[CH1 -] [CH2 -]	<ol style="list-style-type: none"> 1. Reduce output intensity for channel 1. 2. Press to unlock button lock.
[CH1+] [CH2+]	<ol style="list-style-type: none"> 1. Start treatment. 2. Raise output intensity for channel 1 (or 2). 3. First press [CH1-] or [CH 2-] to open button lock (if it has been activated).

 **Caution:** Read the entire cautionary text before use!

7. Placing the electrodes on the skin

1. Make sure that is skin is completely dry and clean.
2. Press the electrodes tightly against the skin so that there is full contact over the entire surface.

Caution!

- Wash, degrease, and dry skin before attaching self-adhering electrodes.
- Do not start the TENS/EMS device before electrodes are firmly in place on the skin.
- Never remove electrodes from the skin when the TENS/EMS device is switched on.

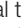
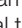
8. Switching on TENS/EMS device


Before you switch on the device we recommend that you read through the cautionary texts, especially under what circumstances the device should not be used. The TENS/EMS device is not a toy but a powerful auxiliary.

Switching on the device: Press and hold down  until the display lights up.

9. Selecting type of treatment

Normal position

- Press  to open normal position. The Display shows:
 - “TENS” – shows that TENS-treatment has been chosen.
 - “P01” – shows that programme 1 (of 20) has been installed.
 - “Clock” – the clock blinks during treatment.
 - “30” – shows the time of treatment selected.
 - “U2” – shows that the treatment applies to user 2.
 - “CH1 0” – shows that channel 1 is not activated (intensity 0). When the channel is activated the indicator lamp CH1 in upper left corner lights up.
- Use pre-installed programme:
Press [] or [] several times to select desired type of treatment. The different types of treatment all have pre-installed programmes.

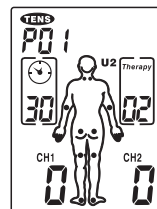EMS (Electrical Muscle Stimulation)	20 different pre-installed programmes for treatment on different places on the body. The different programmes are shown on the display.
TENS (Transcutaneous Electrical Nerve Stimulation)	20 different pre-installed programmes for treatment on different places on the body. The different programmes are shown on the display.
MASSAGE (electric stimulation of muscles in the treatment area)	10 different pre-installed programmes for treatment on different places on the body. The different programmes are shown on the display.

- Activate the programme selected by pressing [CH1+] or [CH2+]. Press several times to increase intensity (0–60).

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Post Clas Ohlson AB, 793 85 INSJÖN
Osoite Clas Ohlson Oy, Annankatu 34–36 A, 00100 HELSINKI