Bathroom scale

Art.no. 44-3815 Model EB1610H

Please read the entire instruction manual before using the product and then save it for future reference. We reserve the right for any errors in text or images and any necessary changes made to technical data. If you have any questions concerning technical problems please contact our Customer Services.

Safety

- This product is a precision instrument and should be handled with care and protected from knocks and rough treatment.
- Place the product in the same place every time you use it.
- Always place the product on an even and level surface for best results. An uneven surface under the scale can cause the readings to fluctuate.
- Weigh yourself without clothes and shoes at the same time every day, preferably in the morning before breakfast. Place your feet in the same place every time you use the scale.
- The displayed readings are rounded up/down to the nearest 100 g.

Batteries


1. Open the battery compartment on the bottom of the scale and insert 2 x AAA/LR03 batteries (sold separately).
2. Refit the battery cover.

Instructions for use


1. Place the scale on an even, hard surface.


2. Stand on the scale and it will come on automatically. The display will stop flashing once weighing is complete.

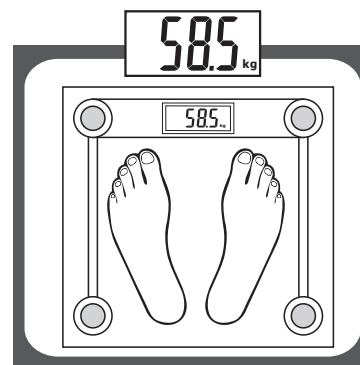3. The scale automatically shuts off as soon as you step off it.

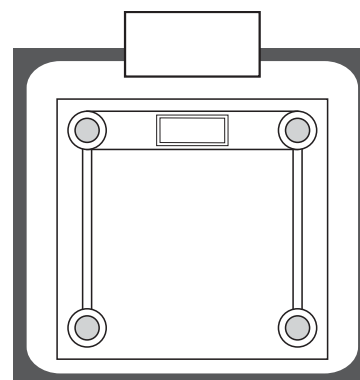4. Press the [Unit button] on the bottom of the scale if you want to change the unit of measurement. You can weigh yourself in kilos (kg) or pounds (lb) or stones (st:lb).

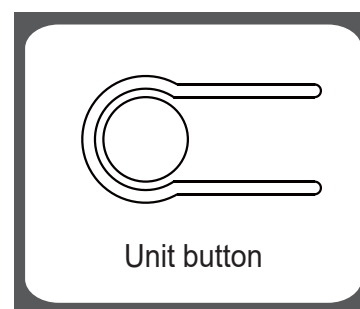Care and maintenance

- Clean the scale by wiping it with a damp cloth. Use only mild cleaning agents, never solvents or corrosive chemicals.
- Never expose the product to liquids or moisture.
- Remove the batteries from the product if it is not to be used for a long time.

The following error messages can appear on the display

Err	The weight capacity of the scale has been exceeded.
Lo	The batteries are low. Insert new batteries.
C	If a C appears on the display when you step off the scale, then weighing was not successful. Weigh yourself again.

If the display does not return to normal or the scale shows an error message, take out the batteries and wait for 5 seconds for the electronics to be reset. Insert the batteries again.

Responsible disposal

This symbol indicates that the product should not be disposed of with general household waste. This applies throughout the entire EU. In order to prevent any harm to the environment or health hazards caused by incorrect waste disposal, the product must be handed in for recycling so that the material can be disposed of in a responsible manner. When recycling your product, take it to your local collection facility or contact the place of purchase. They will ensure that the product is disposed of in an environmentally sound manner.

Specifications

Batteries	2 x AAA/LR03 (sold separately)
Max weight	150 kg
Size	300 x 300 x 20 mm

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